# Smart Body Fat Scale

**SEC 103.02.04**

**SYWB-S007**

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EN – USER MANUAL

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ES – MANUAL DE INSTRUCCIONES

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The screenshots show the following steps:

- Login/Register Screen:** Features a woman in a yoga pose on the left. Buttons for "Register" (green), "Login" (grey), and "Test" (grey) are at the bottom.
- Fitdays Login Screen:** Requests an email address (example@gmail.com). Buttons for "Enter password" (grey) and "Login" (green) are present.
- Enter password Screen:** Requests a password (6 to 16 characters with letters and numbers). Buttons for "Enter password" (grey), "Remember password" (checkbox), and "I forgot my password" (link) are shown.
- Privacy Policy Screen:** Displays the privacy policy text. A checkbox for "I have read and agree to the privacy policy." is checked, and a "NEXT" button is at the bottom.
- Email Registration Screen:** Requests an email address (example@gmail.com). Buttons for "Email Registration" (grey), "Password" (grey), and "Help" (grey) are shown.
- Password Screen:** Requests a password (6 to 16 characters with letters and numbers). Buttons for "Enter password" (grey), "Re-enter password" (grey), and "Re-enter password" (grey) are shown.
- Next Step Screen:** Shows "I already have an account." and a "NEXT" button.

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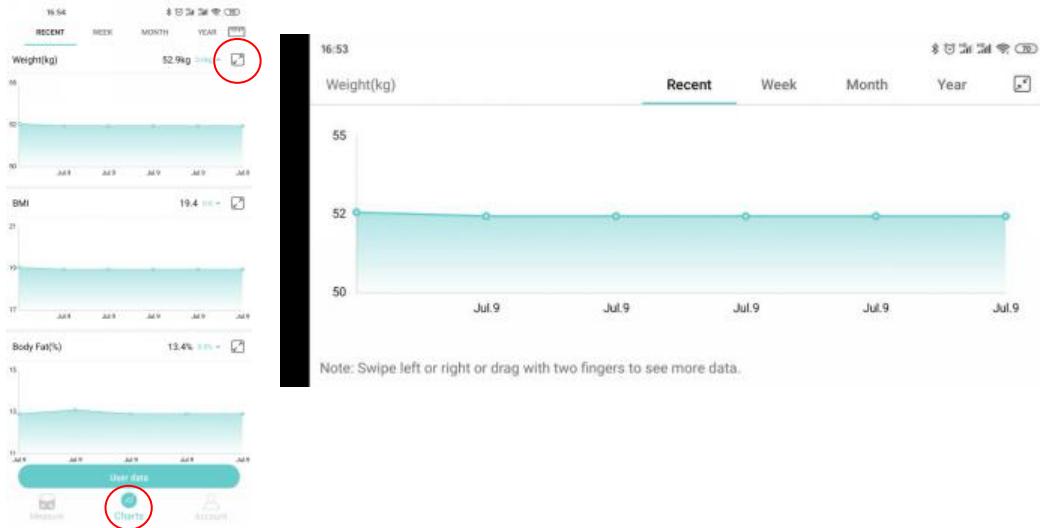
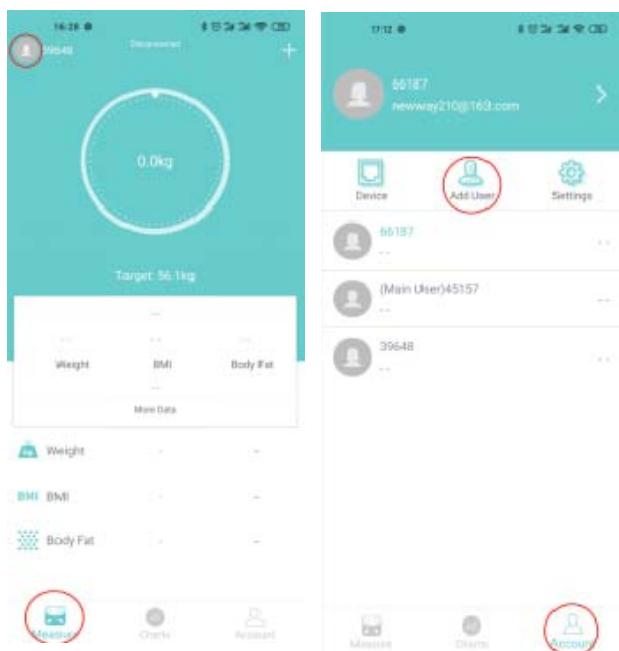
# EN – USER MANUAL

Read carefully this user manual before first use. Pay special attention to the security instructions. Keep this manual for later consultation.

## WARRANTY AND AFTER-SALES

Everyone at SENYA thanks you for your order and for the confidence you have in the Senya brand.

### **Warranty conditions**

As of 1 January 2016, SENYA guarantees the proper functioning of its products for a period of 2 years from the date of purchase. The warranty of any purchase made before that date is valid for 1 year. This warranty only applies to products used for domestic purposes. Any other use in the context of professional use by a natural or legal person excludes any application of the warranty. This product is guaranteed against any failure resulting from faulty workmanship or materials. This warranty does not cover defects or damage resulting from improper installation, incorrect use or normal wear and tear of the product.

### **Availability of spare parts**

Senya has spare parts for its products for a 2-year period. You can find accessories and spare parts for your Senya products sold by our partner Cdiscount. If your accessory or spare part is not available, please contact our after-sales service at [service-client@senya.fr](mailto:service-client@senya.fr).

### **Customer reviews**

If you are satisfied with your product, please submit a review on the website where you placed your order to help other users make their choice. If you are not totally satisfied with your product, please write to us at this email address [service-client@senya.fr](mailto:service-client@senya.fr) before leaving a negative review. Our customer service will listen to your complaints or recommendations in order to satisfy you and improve our products.

### **Bisphenol A Policy**

All our Senya products are developed with great care to meet all European legislation. Our partner laboratories strictly control each production to ensure that our products comply with the requirements for electrical safety, electromagnetic compatibility, limitation of use of certain hazardous substances and food contact for food-related products. All our products are therefore compliant with the current regulations on Bisphenol A (BPA).

## I. SAFETY INSTRUCTIONS

1. Inspect the appliance for visible signs of damage before use. Do not use a damaged device.
2. This appliance is intended for private and non-commercial use. Any use other than previously stated is considered as improper use. No claims of any kind will be accepted for damage resulting from use of the appliance for other than its intended purpose. The risk must be borne solely by the user.
3. Place the product on a stable surface.
4. Under no circumstances must the appliance be placed in water or any liquid or come in contact with liquid.
5. Do not expose the device to extreme temperatures, wide temperature variations, moisture, wet conditions or direct sunlight.
6. Do not drop the device or allow it to be knocked violently. Make sure that no objects can fall onto the scales as the glass could be damaged.
7. This device is intended to be used for determining body weight and analyzing the proportion of water and body fat as well as the proportion of muscle and bone mass in the body, and for calculating the body mass index (BMI).
8. The analysis of the body is not suitable for people under the age of 10 or over the age of 70. We recommend that these people only use the conventional weighing mode.
9. Please note that for technical reasons measurement tolerances are possible as these scales are not calibrated scales for professional or medical use. The results which are calculated are only reference values and act as a guideline. Moreover, different values may occur in the

case of performance athletes or people who are being treated with medicines. You are explicitly advised against using the results as the basis for any form of medical treatment or diet! You must consult your doctor if you have any questions about your weight, a medical treatment or a diet.

10. The set of body analysis scales is not suitable for the following people:
  - People with a pacemaker or another medical implant with electronic components. The correct functioning of the implant could be impaired or endangered.
  - People with symptoms of fever, oedema or osteoporosis
  - People undergoing dialysis treatment
  - People who are taking cardio-vascular medication
  - Pregnant women
11. The device can withstand a maximum body weight of 180 kg ! A higher weight may damage the device!
12. Do not stand on the scales with wet feet and make sure that the measuring surface and the surface underneath is dry. You could slip!
13. Do not stand on the edge of the scales. You could slip off!
14. Keep the device away from electromagnetic fields (mobile phones, remote controls etc.). The measurement accuracy could be impaired.
15. Always place the device on a dry, even and firm surface. Soft surfaces, such as long pile carpets, may distort the measurement result.
16. The device should only be used in closed indoor rooms with a temperate climate. When the device is brought


from a cold room into a warm room, condensation can form. Allow the device to warm back up to room temperature before you use it.

17. Never jump onto the device so that you do not damage the device.
18. Do not use the product in a hot, humid or polluted environment in order to avoid any danger from battery leakage or a short circuit.
19. Avoid any shock or drop of the product.
20. The product is powered by 2 x AAA . Use only fresh batteries of the recommended type.
21. Batteries are to be inserted with the correct polarity (+) and (-) as indicated.
22. Replace all the batteries at the same time.
23. Take out batteries if not use for a long time.
24. Do not dispose of batteries in fire. Batteries might explode or leak.
25. Different types of batteries or new and used batteries are not to be mixed
26. It is advised to supervise the kids that they will not play with the product
27. If the product is damaged, it must be repaired by the manufacturer, its after-sales service, or by a third person with similar qualification, in order to avoid all risks.
28. This appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards

- involved. Children shall not play with the appliance.
29. Waste batteries should not disposed of with household waste. Please check with your local authority for battery recycling advice.

## II. CHARACTERISTICS OF THE PRODUCT

### II.1. GLOBAL DESCRIPTION OF THE PRODUCT
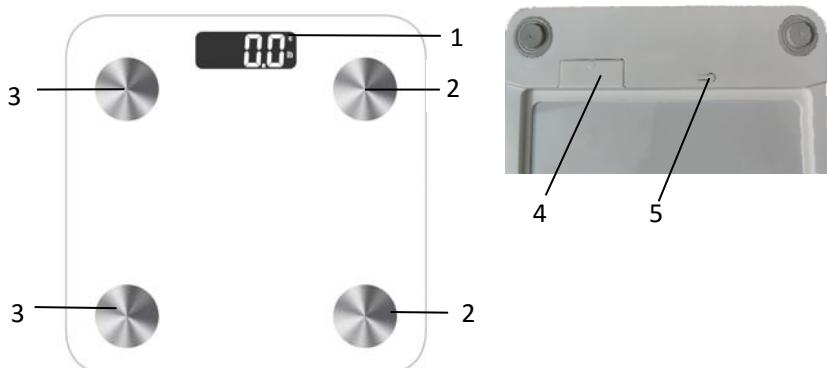
<b>1</b>	LCD Display	<b>2</b>	Contact area for the right foot
<b>3</b>	Contact area for the left foot	<b>4</b>	Battery compartment
<b>5</b>	Unit button		

### II.2. CONTENTS OF PACKAGING

- Smart body fat scale
- User manual

### II.3. TECHNICAL SPECIFICATIONS OF THE PRODUCT

- The scale has the following diagnostic functions: Weight, BMI (Body Mass Index), BFR (Body Fat Rate), Muscle, Moisture, Bone mass, BMR (Basal Metabolic Rate), Visceral fat, Protein rate and Body age and etc.
- Compatible with Android 6.0 or above and iOS 8.0 or above
- Possibility to record up to 24 user profiles
- Weight limit: 6kg-180kg
- Unit: Kg / Lb / St
- Graduation: 100 g
- LCD display: 74x28mm
- Baby scale mode
- Auto power off
- Signal of low voltage (LCD display : LO)
- Powered by 2 x AAA batteries (not included)

- Product size: 260x260x20mm
- Operating temperature: 10-40°C
- Operation humidity: 20%-90%

## **III. USE OF THE PRODUCT**

### **III.1. BEFORE THE FIRST USE**

- Remove the appliance from the packaging and make sure to remove the wedging elements.
- Wipe the external surfaces of the appliance with a soft, clean cloth.

### **III.2. Operation: Conventional weighing**


1. Take out the scale and open up the battery cover, insert the 2 x AAA batteries (not included) into the battery compartment
2. The scale will automatically switch on and will firstly display «88:88» during its initialization. Then the display goes out.
3. Put the scale on a flat steady surface
4. Stand with bare feet on the scales. Make sure that your feet are placed correctly on the contact areas (2 & 3). The body weight will be calculated.
5. Once the body weight has been established, the number on the display begins to flash. You can now get off the scales.
6. The scale will turn off automatically after 10s without any use.

### **III.3. Operation: Connected weighing**

1. Download the app "Fitdays" either using the QR code or from the Google Play Store (Android) or App Store (iOS).
2. Install the App on your device. The app requires iOS 8.0 or higher, or Android 6.0 or higher.
3. The App icon will appear on your phone or tablet after the installation is completed.
4. Switch on the Bluetooth function on your smart phone.


## 5. Register your own ID and Login by email or social media account.

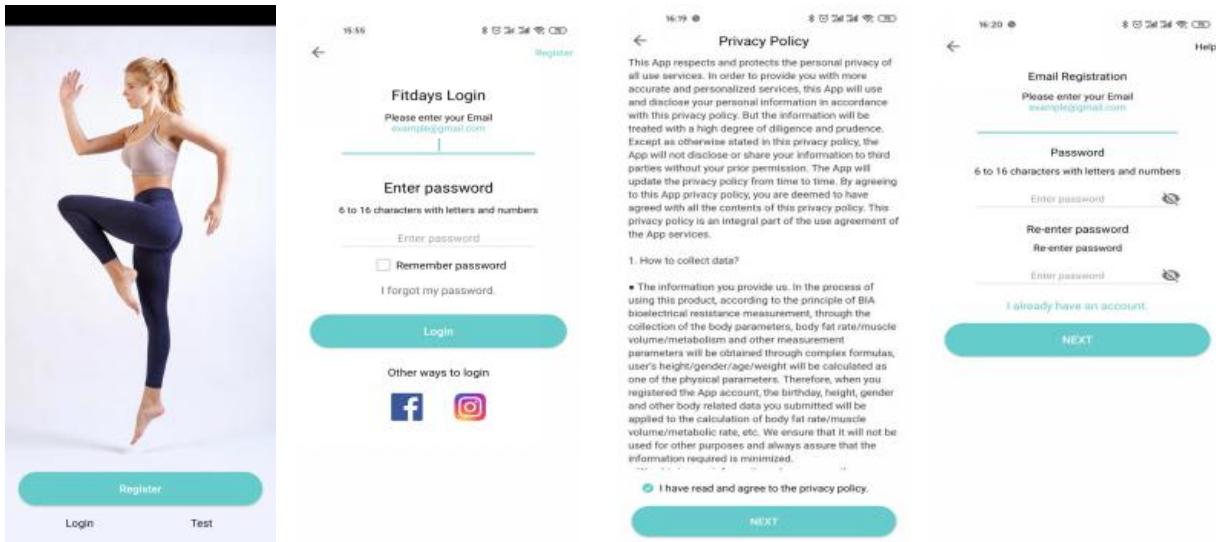
## 6. Add personal data and confirm .

Note: Enter the wrong gender, age, height, or target weight may result in measurement being inaccurate.

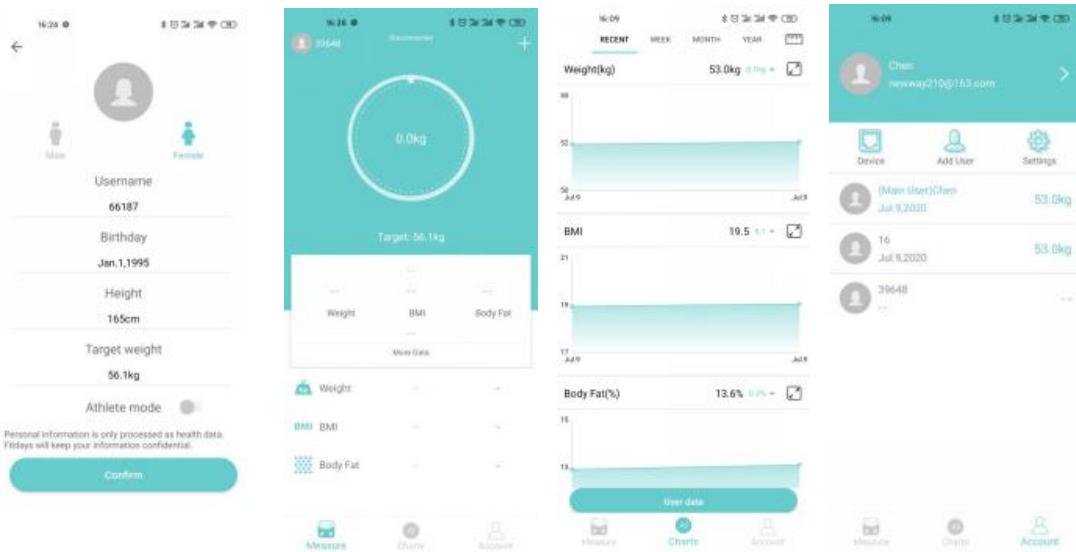
## 7. Pair the scale through Bluetooth. Click: Account ( ) ---> Device ( )---> + ---> find the scale ID pair ( ), now your phone is connected with the scale.

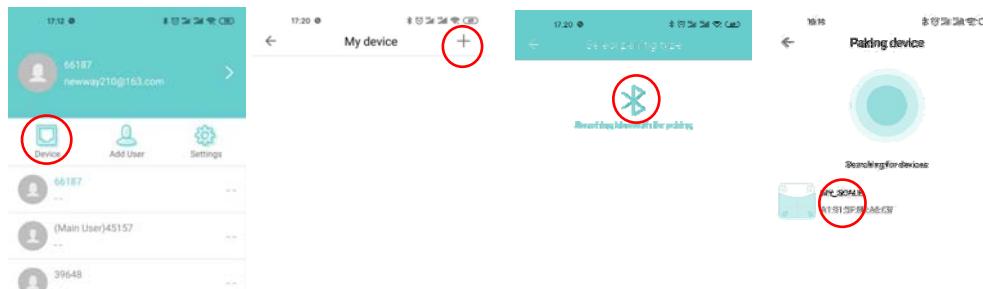
- Tap your foot on the scales to activate the display. After about 2 seconds the scale will show zero "0.0" on the display.
- Stand with bare feet on the scales. Make sure that your feet are placed correctly on the contact areas (2 & 3). The body weight will be calculated.


- The number on the display is flashing when the scale is weighing and analyzing. When LCD show to your weight, You will have about 5-10s to adjust your body ;Continue to stand on the scale until the circle is stable and shows '**Finished**', your phone has received all the data from the scale,You can click "**More Data**" to view various data.

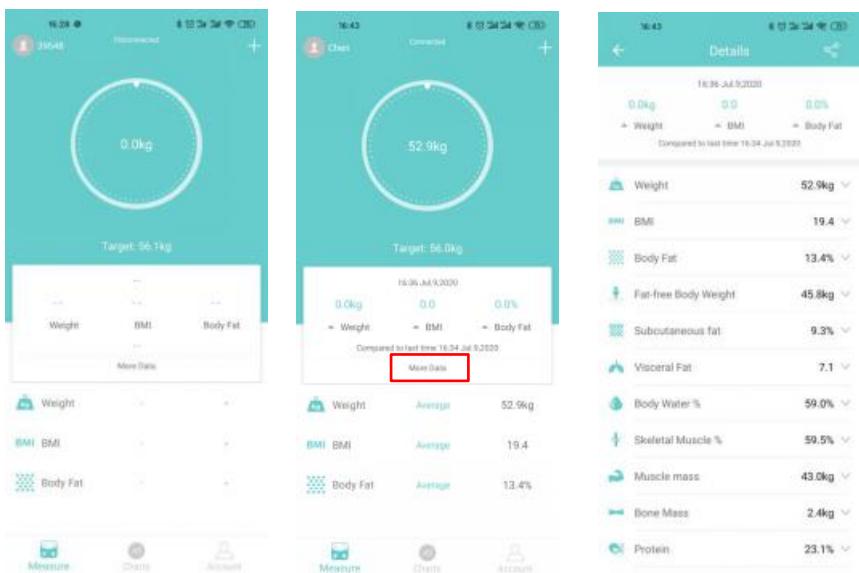- You can also view your daily, weekly, monthly and annual statistics by clicking "**Charts**" and symbol ↗ .

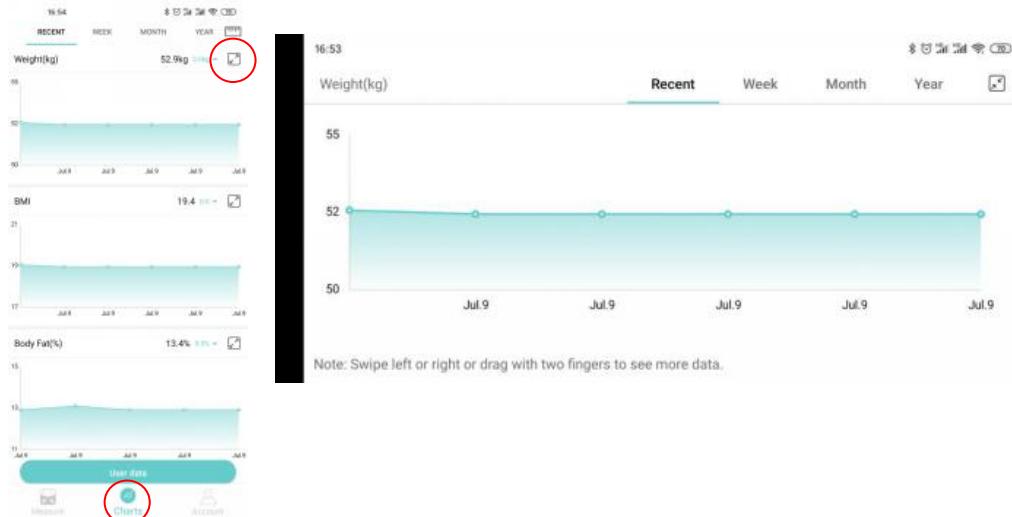12. Add or switch user by the symbol on the upper left (👤). Click **Account** (👤) --> **Add user** (👤) add personal data. The smart scale can recognize which user is step on it and remind user to switch the correct user. Each account can support 24users.

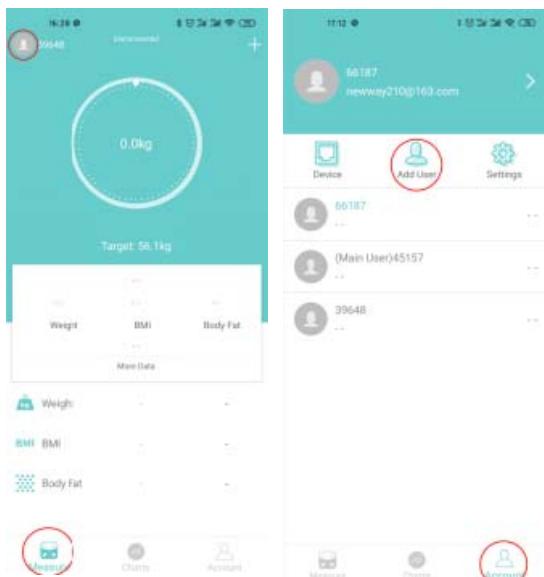13. Weigh baby: click + --> **Baby Mode**. Weigh your weight first— get off the scale — hold the baby and weigh together, the scale will automatically calculate the baby's weight, and finally you can create your baby's account.

## IV.CARE AND STORAGE

1. Put the scale on a clean and flat table.
2. The maximum weight of the scale is 180 kg .
3. To ensure the life of the scale, do not put the scale under wet, hot & cold condition.
4. The scale is a high precision electronic weighing product, please do not press or drop the scale otherwise it is easy to be damaged.
5. If the scale surface is dirty, please clean the scale with soft cloth, do not use any chemical detergent and allow water to get inside the scale.
6. If the result appears improperly, please check whether the batteries and the touching point contact well. Replace for new batteries if used for a long time already, and then try again. If necessary, contact your local distributor or the manufacturer.
7. Do not store the unit near heated sources (eg. Fire, radiators, fan heaters or other heat-emitting devices).

## V. TROUBLESHOOTING GUIDE

1. When the screen shows " Lo " , it means that the power is low, please change the batteries as soon as possible.
2. When the scale is overloaded, " Err " will appear on the display.
3. When an error has occurred during measurement, " C " will appear on the display. Step off the scale and step onto them again to repeat the measuring process.
4. Find the **Unit button** at the back of the scale, press it to exchange unit from Kg to Lb or St.

## V. TREATMENT OF THE ELECTRIC AND ELECTRONIC APPLIANCES AT THE END OF LIFE

Application in the countries belonged to the European Union and in other European countries which have a selective collect system. This symbol on the product or on its packaging indicates that this product shall not be treated as household waste. Instead it shall be handed over to the applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The recycling of materials will help to conserve natural resources. For more detailed information about recycling of this product, please contact your local city office or your household waste disposal service. The packaging material is recyclable. Eliminate the packaging in an environmentally friendly way and make it available to the recycling collection service.

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The image shows a sequence of screenshots from the Fitdays mobile application:

- Login Screen:** Displays a woman in a dynamic yoga pose on the left. The center has fields for "Email" (example@gmail.com) and "Enter password". Below are "Remember password" and "I forgot my password" checkboxes, and a teal "Login" button.
- Privacy Policy Screen:** Shows the "Privacy Policy" section with a detailed text about data collection and processing. It includes a "NEXT" button at the bottom.
- Email Registration Screen:** Requests an "Email" (example@gmail.com) and a "Password" (6 to 16 characters with letters and numbers). It also asks for "Re-enter password" and "I already have an account?" checkboxes. A teal "NEXT" button is at the bottom.

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The first screenshot shows a woman in a dynamic pose on the left, and a 'Fitdays Login' form on the right. It includes fields for email ('Please enter your Email example@gmail.com'), password ('Enter password'), and a 'Remember password' checkbox. Below these are 'Login' and 'Other ways to login' buttons with Facebook and Instagram icons. At the bottom are 'Register' and 'Test' buttons.

The second screenshot shows the 'Privacy Policy' screen with a large amount of text about data protection and user rights. It includes a 'NEXT' button at the bottom.

The third screenshot shows the 'Email Registration' screen with fields for 'Email Registration' ('Please enter your Email example@gmail.com') and 'Password' ('6 to 16 characters with letters and numbers'). It also includes 'Enter password', 'Re-enter password', and 'I already have an account?' links, along with a 'NEXT' button.

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The first screenshot shows a profile setup screen with fields for 'Username' (66187), 'Birthday' (Jan.1,1995), 'Height' (165cm), 'Target weight' (56.1kg), and 'Athlete mode' (radio button). It includes a note about data processing and a 'Confirm' button.

The second screenshot shows a main dashboard with a large circle for weight (0.0kg) and BMI (19.5). It also displays a bar chart for body fat percentage (13.6%). Below these are buttons for 'Weight', 'BMI', and 'Body Fat'.

The third screenshot shows a 'User data' screen with a list of recent measurements for a user named 'Chen'. It includes columns for 'Date', 'Weight' (e.g., 53.0kg), and 'BMI' (e.g., 19.5). At the bottom are navigation buttons for 'Measure', 'Charts', and 'Account'.

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